

NaviCross[®]

Support Catheter

SUPPORTING
YOUR SUCCESS



from ACCESS
to **CLOSURE**
PERIPHERAL
INTERVENTION



**TERUMO**
INTERVENTIONAL
SYSTEMS

Not all products are available for sale in all countries. This information is provided only in respect to markets where these products are approved or cleared. All products are not cleared or approved in the U.S.A. by the Food and Drug Administration. Please contact your Terumo local sales representative for more information.

NaviCross®

Support Catheter

SUPPORTING YOUR SUCCESS

NaviCross® Support Catheter is dedicated to reach and cross complex lesions in challenging anatomy.

Available in 0.035" and 0.018" wire compatibility, NaviCross supports your success in both femoro-popliteal and below-the-knee interventions.

Stainless Steel Double-Braiding

Highly kink-resistant with great torque control and optimal push transition from proximal shaft to distal tip.

Small Crossing Profile

Provides seamless guidewire-to-catheter transition leading to high support and successful lesion crossing.

Hydrophilic M Coat™ 40 cm Distally

Ensures Best-in-Class trackability in challenging anatomy and great crossability of complex lesions.

Three Radiopaque Markers

Facilitates accurate assessment of position and intraluminal measurement of stent and balloon sizes.

NaviCross®

Support Catheter

SUPPORTING YOUR SUCCESS

TWO DEVICES WORKING AS ONE

NaviCross® 0.018" can be smoothly advanced through NaviCross® 0.035" while providing a seamless catheter-to-catheter transition.

This combination offers you:

- More pushability* when tackling difficult stenoses.
- Ability to reach more distal lesions.

RADIFOCUS® GLIDEWIRE ADVANTAGE™ – THE PERFECT PARTNER

NaviCross® and Radifocus® Glidewire Advantage™, the optimal combination for your lower limb interventions by ensuring great support for lesion crossing.

*Compared to the use of 1 device. Data on file.

NaviCross®

Support Catheter

PRODUCT SPECIFICATIONS

Braiding: Double-braided stainless steel

Wire compatibility: 0.035", 0.018"


Outer diameter: 4 Fr (NaviCross® 0.035")
2.6 Fr (NaviCross® 0.018")

Lengths: 65, 90, 135, 150 cm

3 Radiopaque markers: 1 mm from tip - 40 mm and 60 mm from precedent marker


Tip: Straight and 30° Angled

NaviCross® 0.035"

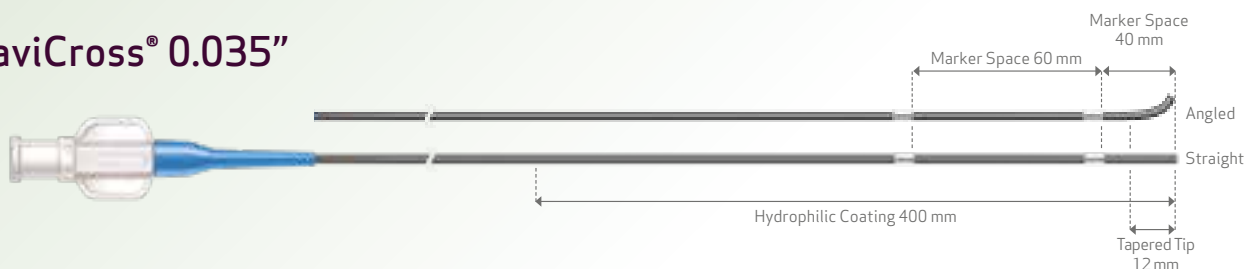

	0.035" / 0.89 mm Wire compatibility	65 cm	90 cm	135 cm	150 cm
Tip shape	Straight	WS*NS350G3HM	WS*NS35093HM	WS*NS350N3HM	WS*NS35153HM
	Angled	WS*NA350G3HM	WS*NA35093HM	WS*NA350N3HM	WS*NA35153HM

NaviCross® 0.018"

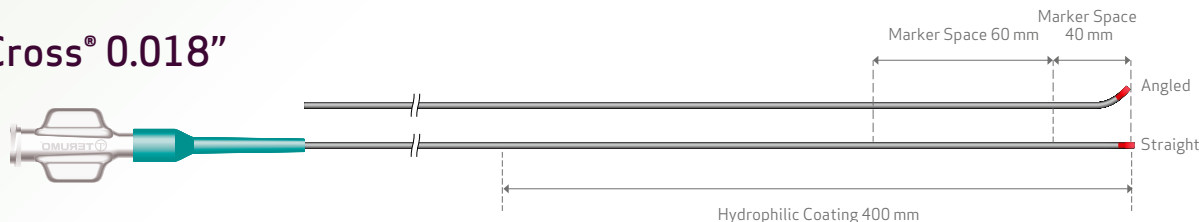

	0.018" / 0.46 mm Wire compatibility	65 cm	90 cm	135 cm	150 cm
Tip shape	Straight	PN*NS180G3HM	PN*NS18093HM	PN*NS180N3HM	PN*NS18153HM
	Angled	PN*NA180G3HM	PN*NA18093HM	PN*NA180N3HM	PN*NA18153HM



Terumo Corporation
+81333748111

Terumo Europe NV
+3216381211

**Terumo
Interventional Systems EMEA**
+33147160930

EMEA Sales Offices

**Terumo Europe NV
Africa Business Division**
+3216381308

**Terumo Europe NV
Benelux Sales Division**
Belgium:
0800 14468
The Netherlands:
0800 0231938

**Terumo Europe NV
Emerging Market Division**
+3216381211

Terumo Deutschland GmbH
+49 6196 80 230

**Terumo Deutschland GmbH
Österreich**
+43 2236 378020

**Terumo Deutschland GmbH
Switzerland**
+41 56 419 10 10

Terumo Europe España SL
+34 902 10 12 98

Terumo France S.A.S.
+33 130 96 13 00

Terumo Italia S.r.l.
+39 06 94 80 28 00

Terumo Russia LLC
+7 495 988 4740

Terumo Sweden AB
+46 3174 85 880

Terumo Middle East FZE
+971 4 292 0200

Terumo UK Ltd
+44 1276 480 440

**Terumo BCT Tıbbi Cihazlar Dağıtım
ve Hizmetleri A.Ş.**
+90 216 645 92 00



© Registered Trademark
Published by Terumo Europe NV

www.terumo-europe.com

Not all products are available for sale in all countries. This information is provided only in respect to markets where these products are approved or cleared. All products are not cleared or approved in the U.S.A. by the Food and Drug Administration. Please contact your Terumo local sales representative for more information.

157926B06181T(0119)E